


# SBBC203 Baby Scale

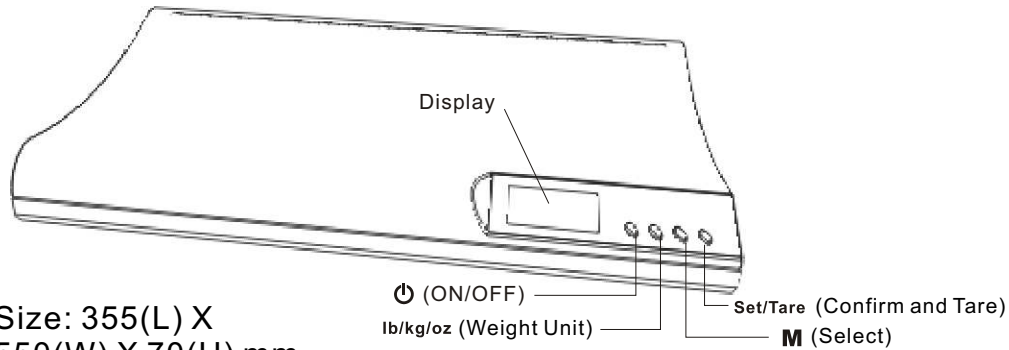

Size: 355(L) X  
550(W) X 70(H) mm


## Features:

- Maximum weight capacity: 20kgs/ 44lbs/ 704oz
- Graduation: 0.01kg/ 0.02lb/ 0.5oz
- Color digital readout
- Memory position: 5 memories in each cartoon character
- Low battery indicator
- Power: 4 X AAA batteries


## Preparations:

1. Battery Installation: open the battery cover at the bottom of the scale and insert 4 brand-new AAA batteries properly.
2. Operating Environment: lay the scale onto a firm and horizontal ground.

## Weighing a naked baby

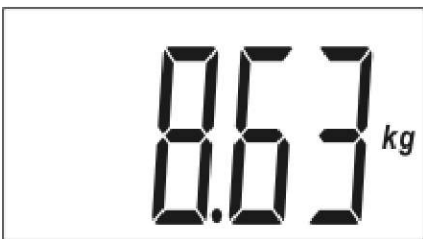


1. Press ON/OFF button to switch on.
2. Last record shows with weight unit in kg. To review the weight unit in lb or oz, press lb/kg/oz button.

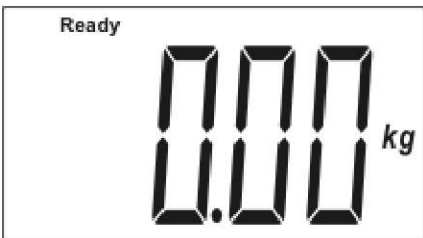

3. After a few seconds, 0.00 and weight unit shows. "Ready" icon appears.  
\*Press lb/ kg/ oz button to select weight unit if different unit is required.

1.

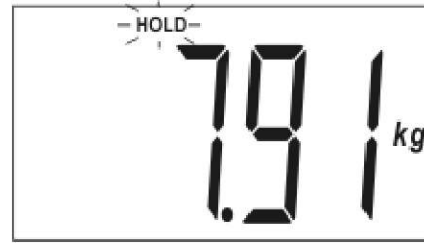


4. Put the baby onto the scale. When "HOLD" flashes, record of the baby is stored.  
\*If a second baby is waiting to be weighed, lift the baby from the scale when "HOLD" is flashing. "----" appears and change to 0.00 (as shown on diagram 2). Put the second baby onto the scale.


## Weighing with towel or container




1. Press ON/OFF button to switch on.
2. Last record shows with weight unit in kg. To review the weight unit in lb or oz, press lb/kg/oz button.

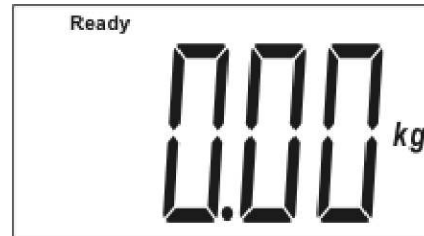

3. After a few seconds, 0.00 and weight unit shows. "Ready" icon appears.  
\*Press lb/ kg/ oz button to select weight unit if different unit is required.


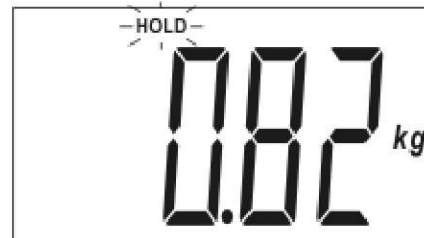


4. Put the towel or container onto the scale.
5. Press Set / Tare button when the result of the towel or container is shown with "HOLD" icon flashing.

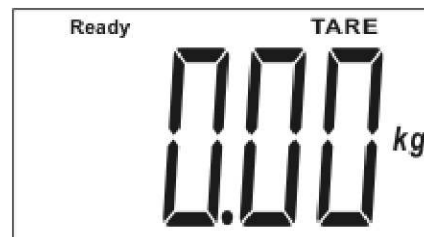


6. Put the baby onto the towel or container when 0.00, weight unit, "Ready" icon and "TARE" icon show.

2.

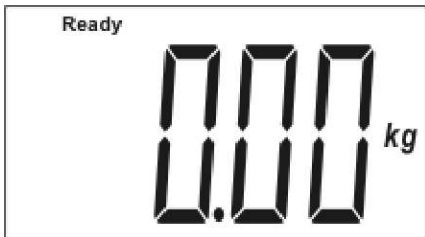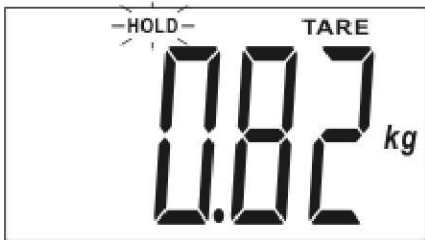

7. Result of the baby will show with "HOLD" icon flashing. After a few seconds, "----" appears and the display will change to standby mode (as shown on the diagram). Power will go off automatically after 50 seconds if there is no second baby put onto it. Power can also be turned off by pressing ON/OFF button.

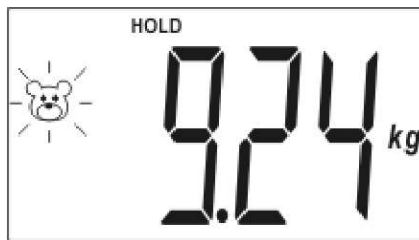

9. Press "Set / Tare" button to confirm. Selected character will flash and record of the baby is stored. Power will go off automatically after a few seconds. When the same baby is weighed next time, repeat steps 1-7. Press "M" button to go to the cartoon character selected before and press "Set / Tare" to confirm. A second record is stored. Power will go off automatically after a few seconds.


## M Button

The main function of the "M" Button is selection and memory. For easy memory, this scale is designed with four different cartoon characters (duck, bear, elephant and dog) to store weight records. Each character can store 5 records. When the sixth record is stored, the first record will be out.


### Weighing with "M" button

Repeat steps 1,2,3 if only a naked baby is weighed.

Repeat steps 1,2,3,4,5,6 if a baby with towel or container is weighed together.

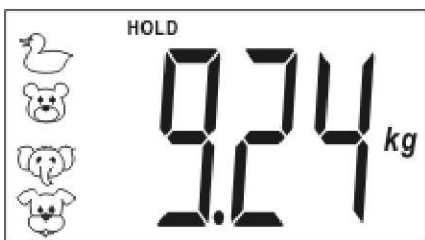


7. When the result of the baby shows with "HOLD" icon flashing, press "M" button. Cartoon characters and result of the baby will display.

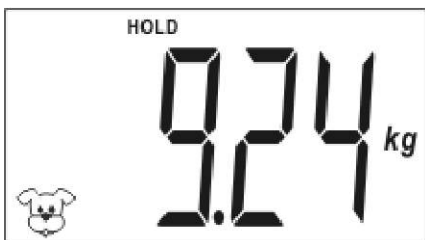


8. Press "M" button again to select cartoon character.

To review the records of the baby stored (only when the power is off):

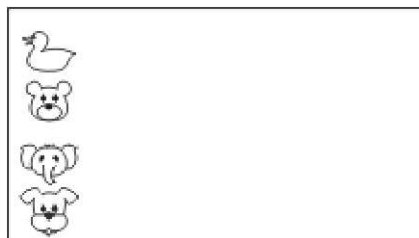


1. Press "M" button to switch on. Four cartoon characters appear.

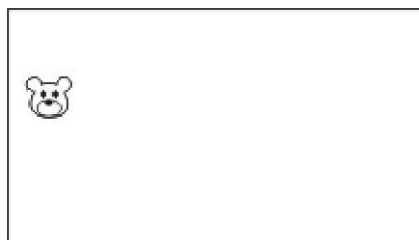

2. Press "M" button to switch to desired character.

3. Press "Set / Tare" button.

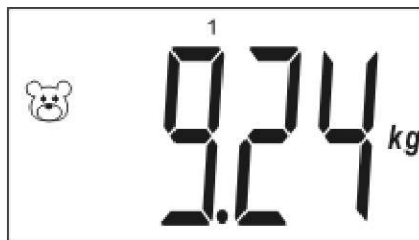

4. Record number 1 will show.

5. Press "M" button again to see different records.

6. Press ON/OFF to switch off. Power will go off automatically in 1 minute if no other button is pressed.



Web site : <http://www.biocareplus.com>  
E-mail : [info@biocareplus.com](mailto:info@biocareplus.com)